

IPL Hair Remover

User Manual




DW DARTWOOD

About the IPL Hair Remover

The Dartwood IPL (Intense Pulsed Light) Hair Remover gets rid of unwanted body hair by directing intense light at the hair root.

In the Box

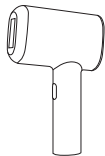


1 x IPL
hair remover



1 x safety glasses



1 x razor

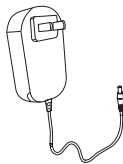


1 x
power adapter

Safety Information

- Read this manual and the instructions carefully before using.
- Failure to follow the safety instructions may result in injury.
- Do not use this device if you are under the age of 16.
- Do not use it if the skin of the area you want to treat is sunburned.
- Do not use it if you are pregnant, trying to become pregnant, or breastfeeding.
- Do not use it if you have epilepsy, get light-induced seizures, or migraines.
- Do not use on tattoos, freckles, birthmarks, moles, or warts.
- Do not use in sensitive areas such as your head, ears, genitals, etc.
- Do not use if you have any active skin conditions in the treatment area.
- Do not use if you have a history of scars.
- Do not use if you have a disease related to light sensitivity.
- Clean the skin treatment area thoroughly.
- Use skin moisturizer or sunblock after treatment.
- Use a cloth or cotton ball to clean treatment window.
- Use a dry cloth or alcohol-free damp cloth to clean the exterior.
- Keep this device away from water and other liquid.
- Store in a ventilated and dry area.
- Do not insert metal objects into the device.
- Do not disassemble or try to repair the device.
- Contact the manufacturer if the device fails to function properly.

Introduction to IPL Hair Remover

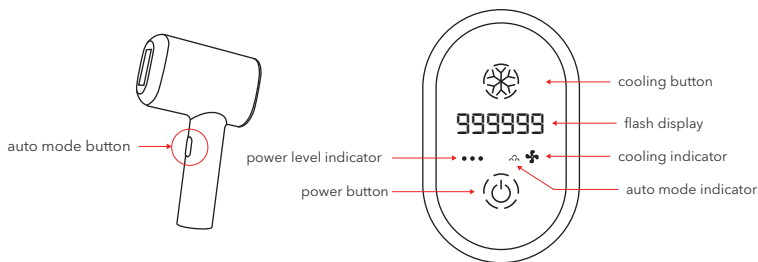

Before You Use

CAUTION:

Always wear the safety glasses while using this device. Avoid looking directly at the light.

SKIN SENSITIVITY TEST:


Test your skin sensitivity by placing the device over a small area of skin at level 1 intensity. Then, move to other areas and gradually increase the intensity.

Wait at least 60 seconds to make sure there are no adverse effects before you continue (redness, irritation, discomfort, etc).


How to Use

1. Clean your skin and shave the hair around the area that you want to treat.


Note: It's best to shave about 2 hours before treatment.

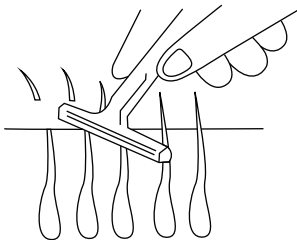


2. Clean the lamp window on the hair remover and plug it in.

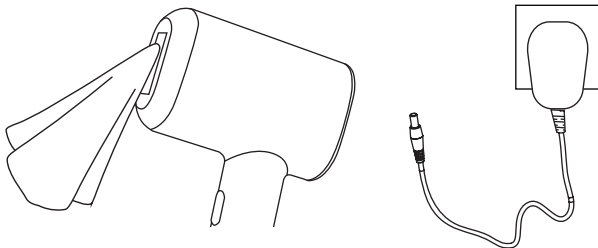


3. Put on the safety glasses to protect your eyes from the IPL light.

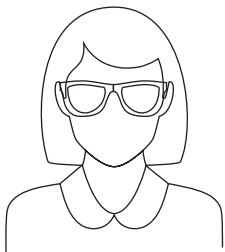

4. Press and hold the power button on the hair remover. Then, lightly press the power button again to adjust the intensity from level 1 - 5.

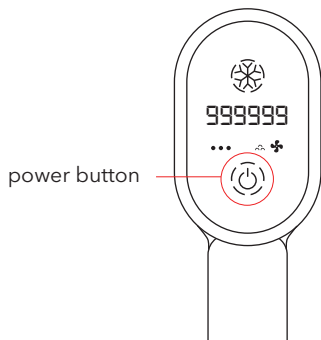


5. Press the cooling button to activate the cooling mode (recommended).

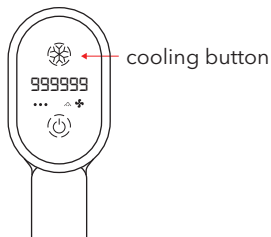


6. Place the hair remover directly on the skin you wish to treat. To activate the pulsing light, press the trigger on the handle (one flash per press). Move slowly around the whole area.

Caution: To avoid burning your skin, do not overlap.

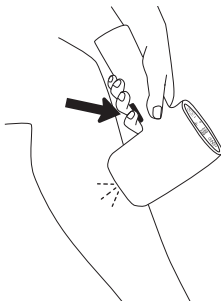

7. Power off and unplug the adapter after use.
8. After treatment, moisturize treated area in 15 minute intervals.

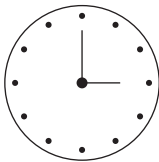

Note: Mild skin redness after treatment is normal. Wait at least two months between treatments in the same area.

Treatment Recommendations

1. For the first 3 treatments, do one treatment every week.
2. For 4 - 9 treatments, do one treatment every two weeks.
3. For maintenance: Do one treatment every two months (in any area with hair regrowth).

Note: Device is most effective at removing black, brown, dark brown, and blonde hair from lighter skin tones. For sensitive or darker skin tones, use a lower energy level.

Possible Side Effects

Common Side Effects	Action
Minor skin irritation or redness in the treatment area	Apply an ice pack or aloe vera lotion to soothe the skin. If skin remains irritated for more than 48 hours, seek medical advice before carrying out further treatment.
Rare Side Effects	Action
Severe redness or swelling	Apply an ice pack or aloe vera lotion to soothe the skin. Seek immediate advice before carrying out any further treatment.
Very Rare Side Effects	Action
Burns or blistering of the skin	Stop treatment and seek immediate medical attention.
A temporary change in skin color	Stop treatment and seek medical attention.

Specifications

Emitted light:	Intense Pulse Light
Emitted energy:	18J max
Energy Level:	5 Levels
	Reference Energy Test Data
	Level 1: 7.4J
	Level 2: 9.2J
	Level 3: 10.7J
	Level 4: 12.5J
	Level 5: 14.4J
Treatment area:	1.7in ²
Input power:	36W
Voltage rating:	100V-240V
Rated frequency range:	50 Hz/60Hz
Storage temperature:	32 °F - 113°F
Working temperature:	41 °F - 95°F
Working humidity:	25% - 75%
Product size:	7.36 x 4.53 x 1.81in
Net weight:	0.52lb



CONTACT US



For more information, please visit:

www.wasserstein-home.com

or contact us at: contact@wasserstein-home.com

Made in China